
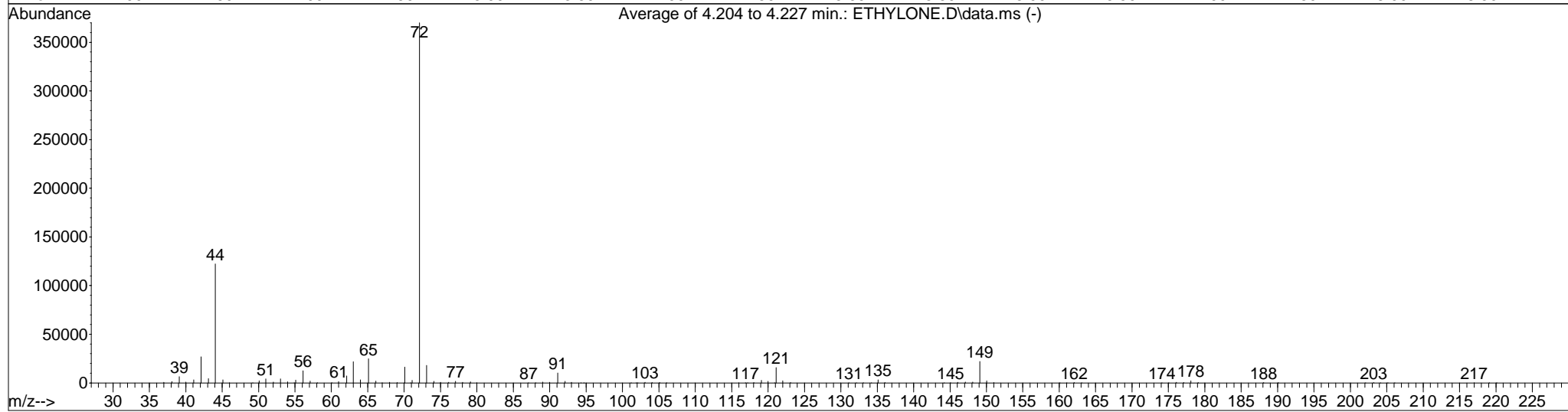


Name of File: D:\Data\forendex\ETHYLONE.D
Method: DRUG1.M
Misc Info: Cayman lot 0432915-10 / DEA SWL

Vial Number: 7
Instrument: DEA# 361782

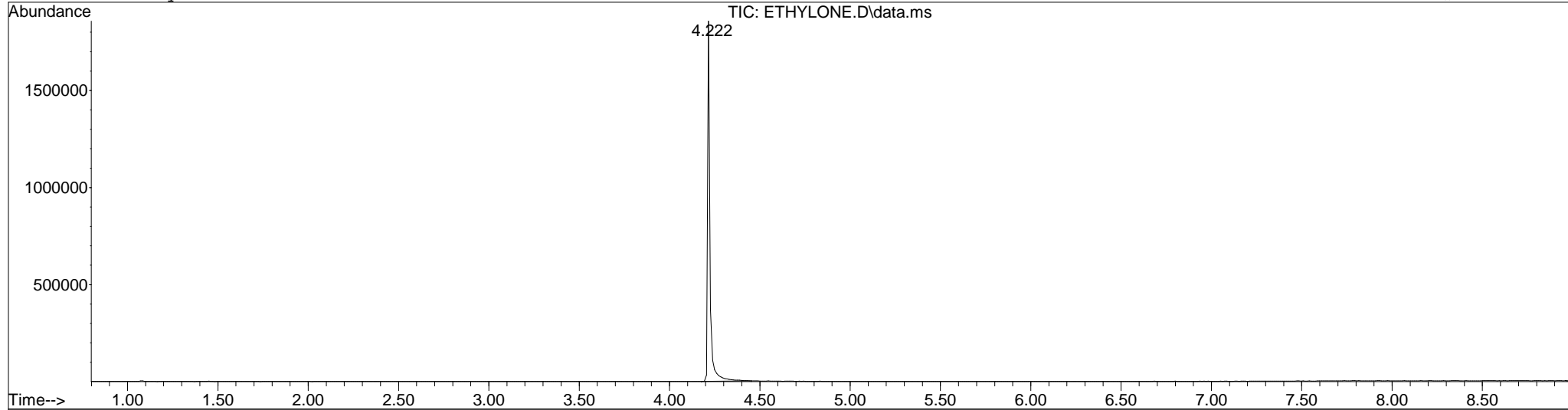

INSTRUMENT CONDITIONS

temp: 90 °C for 0.5 min
then 30 °C/min to 300 °C for 1.5 min
Run Time 9 min
inj. port: 280 °C, transfer line: 280 °C
GC column: ZB-5MS: 15 m x 250 µm x 0.25 µm film
split mode: Split, ratio: 50 :1
tune: stune.u , multiplier offset: 0V

Laboratory Number: Ethylone in meoh
Date: 17 Feb 2012 17:53
Initials: _____